


Ovarian shield 278

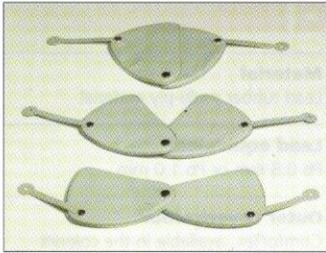

Material

Lead rubber multi-ply material

Lead equivalent

Pb 1.0 mm

Outer cover/colour

Soft PVC pale-grey colour

Sizes


RP2781 – large (12 cm high)
RP2782 – medium (9 cm high)
RP2783 – small (7 cm high)
Available also as complete set (all 3 sizes)
at a lower price = RP278SET!

The design of this MAVIG ovarian shield is unsurpassed and meets the requirements of radiology practice, past and present. The radiation load on the female gonads is reduced considerably without interfering with the diagnostic techniques.

Extensive studies provided the basis for the size and position relation between the anterior iliac vertebrae (fixed point), the ovary shield, and its projection on the pelvic image. The two control marks that are permeable to radiation facilitate proper and uncomplicated positioning.

Each of the three sizes is adjustable thanks to the fan shape that allows for optimum fine tuning. A set with all sizes offers the best possible protection for all of your patients.

It goes without saying that the shield is easy to clean and disinfect.

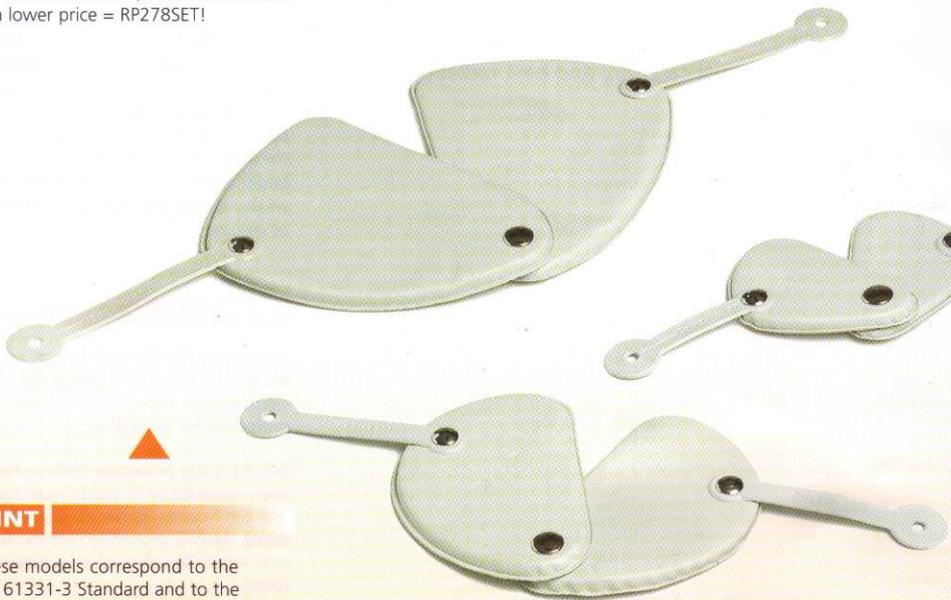

HINT

These models correspond to the IEC 61331-3 Standard and to the German rules concerning Radiation protection for medical use of X-Rays.

MAVIG X-Ray protective clothing has been approved by EU model examination and certification, as per EU Guideline 93/42. MAVIG protective clothing fulfill with regard to safety all legal requirements.